

# Mouse IgG3 ELISA Kit

## Life Diagnostics, Inc., Catalog Number: 5010-1D

### ELISA for the Quantitative Determination of Mouse IgG3 in Serum or Plasma

#### INTRODUCTION

The mouse IgG3 ELISA kit is designed for measurement of IgG3 in mouse serum or plasma. The assay uses goat anti-mouse IgG for solid phase (microtiter wells) immobilization and horseradish peroxidase (HRP) conjugated goat anti-mouse IgG3 antibodies for detection. The capture antibodies were cross-adsorbed on mouse IgM and IgA agarose columns, and the detection antibodies were cross adsorbed on mouse IgM, IgA, IgG1, IgG2a and IgG2b agarose thereby ensuring specificity for IgG3. Cross-reactivity with immunoglobulins from other species has not been investigated.

#### PRINCIPLE OF THE TEST

Test samples are diluted and incubated in the microtiter wells for 45 minutes alongside prepared mouse IgG3 standards. The microtiter wells are subsequently washed and HRP conjugate is added and incubated for 30 minutes. IgG3 molecules are thus sandwiched between the immobilization and detection antibodies. The wells are then washed to remove unbound HRP-labeled antibodies and TMB Reagent is added and incubated for 20 minutes at room temperature. This results in the development of a blue color. Color development is stopped by the addition of Stop Solution, changing the color to yellow, and optical density is measured spectrophotometrically at 450 nm. The concentration of IgG3 is proportional to the optical density of the test sample and is derived from a standard curve.

#### MATERIALS AND COMPONENTS

##### **Materials provided with the kit:**

- Anti mouse IgG coated 96-well plate (12 strips of 8 wells)
- Anti mouse IgG3-HRP Conjugate, 11 ml
- Reference standard (lyophilized)<sup>1</sup>
- 20x Wash Solution, 50 ml
- 10x Immunoglobulin Diluent (25 ml)
- TMB Reagent (One-Step) 11 ml
- Stop Solution (1N HCl), 11 ml

##### **Materials required but not provided:**

- Precision pipettes and tips
- Distilled or deionized water
- Polypropylene or glass tubes
- Vortex mixer
- Absorbent paper or paper towels
- Micro-Plate incubator/shaker mixing speed of ~150 rpm
- Plate washer
- Plate reader with an optical density range of 0-4 at 450nm
- Graph paper (PC graphing software is optional)

<sup>1</sup> The reference standard consists of lyophilized mouse serum of known IgG3 concentration in a BSA matrix. The IgG3 concentration was determined relative to purified mouse IgG3 obtained from an independent laboratory.

#### STORAGE OF THE TEST KIT

The test kit will remain stable for six months from the date of purchase provided that the components are stored as described above. The microtiter plate should be kept in a sealed bag with desiccant to minimize exposure to damp air.

#### GENERAL INSTRUCTIONS

1. Please read and understand the instructions thoroughly before using the kit.
2. All reagents should be allowed to reach room temperature (18-25°C) before use.
3. Optimum results are achieved if, at each step, reagents are pipetted into the wells of the microtiter plate within 5 minutes.

#### IMMUNOGLOBULIN DILUENT PREPARATION

The diluent is provided as a 10x stock. Prior to use estimate the final volume of diluent required for your assay and dilute one (1) volume of the 10x stock with nine (9) volumes of distilled or deionized water.

#### WASH SOLUTION PREPARATION

The wash solution is provided as a 20x stock. Prior to use dilute the contents of the bottle (50 ml) with 950 ml of distilled or deionized water.

#### STANDARD PREPARATION

1. The mouse IgG3 standard is provided as a lyophilized stock. Reconstitute as detailed on the vial label (***the reconstituted standard is stable at 4°C for one week but should be aliquoted and frozen at -20°C after reconstitution if future use is intended***).
2. Label 6 polypropylene or glass tubes as 25, 12.5, 6.25, 3.13, 1.56 and 0.78 ng/ml.
3. Into the tube labeled 25 ng/ml, pipette the volume of diluent detailed on the IgG3 standard vial label. Then add the indicated volume of IgG3 standard (shown on the IgG3 standard vial label) and mix gently. This provides the 25 ng/ml standard.
4. Dispense 250 µl of diluent into the tubes labeled 12.5, 6.25, 3.13, 1.56 and 0.78 ng/ml.
5. Prepare a 12.5 ng/ml standard by diluting and mixing 250 µl of the 25 ng/ml standard with 250 µl of diluent in the tube labeled 12.5 ng/ml.
6. Similarly prepare the 6.25, 3.13, 1.56 and 0.78 ng/ml standards by serial dilution.

#### SAMPLE PREPARATION

**General Note:** Mouse serum IgG3 concentrations vary with species and treatment conditions. In our laboratory we have measured values of 0.05 to 0.5 mg/ml or greater. We suggest an initial dilution factor of 50,000 fold using the following

**procedure for each sample to be tested. We also advise that a limited number of samples be tested initially (as singlets) in order to determine the optimal dilution factor.**

1. Dispense 498  $\mu\text{l}$  and 796  $\mu\text{l}$  of 1x diluent into separate tubes.
  2. Pipette and mix 2  $\mu\text{l}$  of the serum/plasma sample into the tube containing 498  $\mu\text{l}$  of diluent. This provides a 250 fold diluted sample.
  3. Mix 4  $\mu\text{l}$  of the 250 fold diluted sample with the 796  $\mu\text{l}$  of diluent in the second tube. This provides a 50,000 fold dilution of the sample.
  4. Repeat this procedure for each sample to be tested
- Tissue extracts, and tissue culture medium will likely have lower IgG3 levels than those found in serum. Optimal dilutions of such samples should be determined empirically.

### ASSAY PROCEDURE

1. Secure the desired number of coated wells in the holder.
2. Dispense 100  $\mu\text{l}$  of standards and diluted samples into the wells (we recommend that samples be tested in duplicate).
3. Incubate on an orbital micro-plate shaker at 100-150 rpm at room temperature (18-25°C) for 45 minutes.
4. Aspirate the contents of the microtiter wells and wash the wells 5 times with 1x wash solution using a plate washer (400  $\mu\text{l}$ /well). The entire wash procedure should be performed as quickly as possible.
5. Strike the wells sharply onto absorbent paper or paper towels to remove all residual wash buffer.
6. Add 100  $\mu\text{l}$  of enzyme conjugate reagent into each well.
7. Incubate on an orbital micro-plate shaker at 100-150 rpm at room temperature (18-25°C) for 30 minutes.
8. Wash as detailed in 4 to 5 above.
9. Dispense 100  $\mu\text{l}$  of TMB Reagent into each well.
10. Gently mix on an orbital micro-plate shaker at 100-150 rpm at room temperature (18-25°C) for 20 minutes.
11. Stop the reaction by adding 100  $\mu\text{l}$  of Stop Solution to each well.
12. Gently mix. *It is important to make sure that all the blue color changes to yellow.*
13. Read the optical density at 450 nm with a microtiter plate reader within 5 minutes.
- 14.

### CALCULATION OF RESULTS


1. Calculate the average absorbance values ( $A_{450}$ ) for each set of reference standards and samples.
2. Construct a standard curve by plotting the mean absorbance obtained from each reference standard against its concentration in ng/ml on linear graph paper, with absorbance values on the vertical or Y-axis and concentrations on the horizontal or X-axis.
3. Using the mean absorbance value for each sample, determine the corresponding concentration of IgG3 in ng/ml from the standard curve.
4. Multiply the derived concentrations by the dilution factor to determine the actual concentration of IgG3 in the sample.
5. PC graphing software may be used for the above steps.
6. If the  $OD_{450}$  values of samples fall outside the standard curve, samples should be diluted appropriately and re-tested.

### TYPICAL STANDARD CURVE

A typical standard curve with optical density readings at 450 nm on the Y axis against IgG3 concentrations on the X axis is shown below. This curve is for the purpose of illustration only and should not be used to calculate unknowns. Each user should obtain his or her data and standard curve in each experiment.

IgG3 (ng/ml)	Absorbance (450 nm)
25	1.875
12.5	1.280
6.25	0.810
3.13	0.498
1.56	0.326
0.78	0.257

### Typical Mouse IgG3 Standard Curve

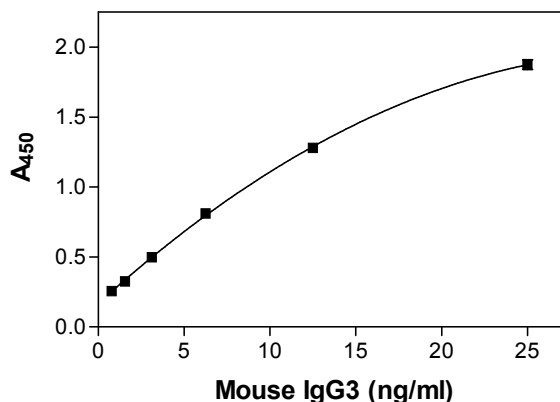

### LIMITATIONS OF THE PROCEDURE

1. Reliable and reproducible results will be obtained when the assay procedure is carried out with a complete understanding of and in accordance with the instructions detailed above.
2. The wash procedure is critical. Insufficient washing will result in poor precision and falsely elevated absorbance readings.